



Brewer
more possibilities

INCANDESCENT SERIES Exam Light

- 1 Crisp white light with no heat
- 2 Advanced reflector system provides shadow-free light
- 3 Flexible, "no-drift", 18" gooseneck offers superior beam alignment
- 4 540° rotation of the light head and gooseneck arm
- 5 Non-metallic head eliminates surface heat
- 6 Internal power cord exits the weighted base at floor level

- 25-Watt incandescent bulb included
- Bulb life: 500 hrs; light pattern: 7" @ 24"
- Maximum height: 62"
- Base diameter: 13"
- Cord length: 9' 3"
- Available with five-caster base
- UL certified & CAN/CSA-C22.2 No.601.1 certified



#15100 - white
#17100 - beige (shown)


accessories

#001525 25 watt replacement bulb
#100165 five-caster base - white
#200165 five-caster base - beige

CHROME SERIES Gooseneck Exam Lamp

- 1 Durable chrome finish
- 2 12 foot power cord w/ grounded 3-wire plug
- 3 Stable chrome plated base (11" diameter)
 - Height: 47" - 72"
 - 10 to 40-Watt incandescent bulb (not included)
 - Some assembly required
 - UL certified

#11500

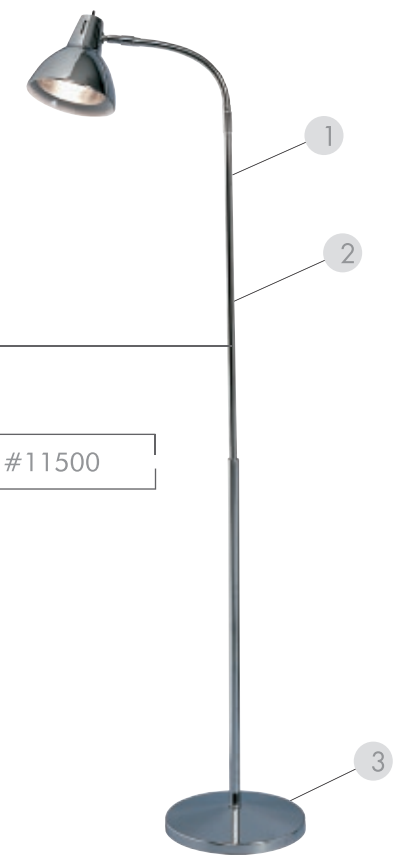

All Brewer lighting comes with a 1-year warranty.
The Brewer Company reserves the right to improve or modify the design of its products consistent with new methods and techniques.

Call your local Brewer sales representative to learn about our complete product line.

The Brewer Company
13901 Main Street, Menomonee Falls, WI 53051
1.888.Brewer.1 www.brewercompany.com

Brewer
more possibilities

FEATHERLIGHT

SERIES Delivers unequal flexibility with one-hand, leading counter-balanced movement.



Dichroic film-coated reflectors filter out 80% of beam heat

Three color-corrected, low voltage, 50-Watt halogen bulbs w/ heat sink technology that keeps the housing cool

"No-drift," counter-balanced arm design; 540° rotation at both the light head and yoke
Arm length: 36.6"

- Ideal for dermatology, plastic surgery, and minor procedures in the clinical environment.
- Bright, crisp light with no heat
- Movement provides effortless and accurate light placement without drifting or repositioning
- 5000 Footcandles @ 24"
- Color temperature - 4000° Kelvin
- Fused for your protection
- Bulb life: 4000 hrs; light pattern: 10" @ 24"
- Floor/Wall/Ceiling models
- Removable and sterilizable handle
- UL certified & CAN/CSA-C22.2 No.601.1 certified
- Ceiling models maximum length: 8': 25"; 9': 37"


Diameter: 25"
Height: 5"
Compact five-leg caster base; power cord exits the base at floor level

accessories


#016150 50-watt replacement bulb



#163002

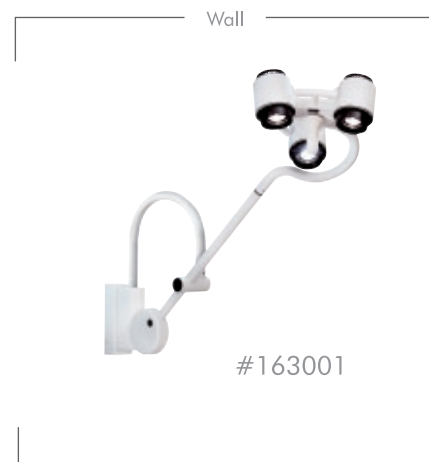

#163001

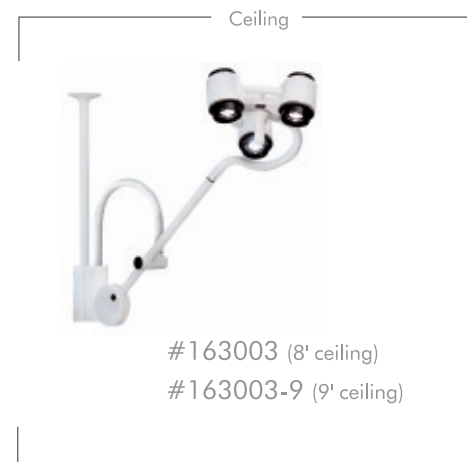

#163003 (8' ceiling)
#163003-9 (9' ceiling)

HALOGEN 35

SERIES Offers flexibility and portability in one superior product.

- Crisp white light with no heat
- Non-metallic head eliminates surface heat
- High intensity, low voltage, 35-Watt bulb included
- 1500 Footcandles @ 24"
- Color temperature - 3000° Kelvin
- Height: 65.5"
- Base diameter: 13"
- Cord length: 9' 3"
- Bulb life: 4000 hrs; light pattern: 7" @ 24"
- UL certified & CAN/CSA-C22.2 No.601.1 certified

Dichroic film-coated reflectors filter out 80% of beam heat

Heat sink technology removes bulb socket heat and keeps the housing cool

540° rotation of the light head and gooseneck arm

Flexible, "no-drift," 24" gooseneck

#18100 beige

accessories

#200165 five-caster base - beige

#001635 35-watt replacement bulb

Internal power cord exits the weighted base at floor level

Shown with optional five-caster base

